

Slicer RegistrationQA Introduction

A registration quality assurance module for 3D Slicer



Universitätsklinikum
Erlangen



Dr. Kristjan Anderle

k.anderle@gsi.de

Introduction

Slicer RegQA extension enables quality assurance (QA) for image registration. It can perform different test, both qualitative and quantitative to estimate registration quality. It can also create output files, which can serve as documentation for specific registration QA.

Terminology

To provide a clear and consistent description of the module, an overview of the nomenclature used is given here.

- *Fixed image* – the image that serves as a reference position in the registration (image that is being registered to).
- *Moving image* – the image that is matched to the reference image (image that is being registered from).
- *Reference image* – the initial choice of fixed image. The choice is based on the registration usage (i.e. opposite registration directions are need for contour propagation and dose wrapping). In a 4D-CT registration a single phase is usually chosen as a reference phase.
- *Transformation* – the result from registration algorithm is a transformation that transforms moving to fixed image. It can either be a matrix for rigid/affine registrations or a vector field for deformable registrations.
- *Warped image* – the result of applying a transformation to the moving image. It should be as close to the fixed image as possible.
- *Forward registration* – registration when reference image is the fixed image. The term “forward” can also be applied to corresponding items (warped image, transformation, Jacobian determinant).
- *Backward registration* – registration when reference image is the moving image. The term “backward” can also be applied to corresponding items.

Measures

The registration QA is based on the recent AAPM task report [1], which recommends verification of 8 different measures. The reader of this document is encouraged to read the document [1], since it provides detailed instructions on registration QA.

Eight different inputs are required to evaluate all measures for a single registration. Measures 1 – 4 are qualitative that require user to spot any irregularities in registration, while measures 5 – 8 are quantitative, where each measure has a tolerance value. All measures, except measure 8, can be verified for both forward and backward registration.

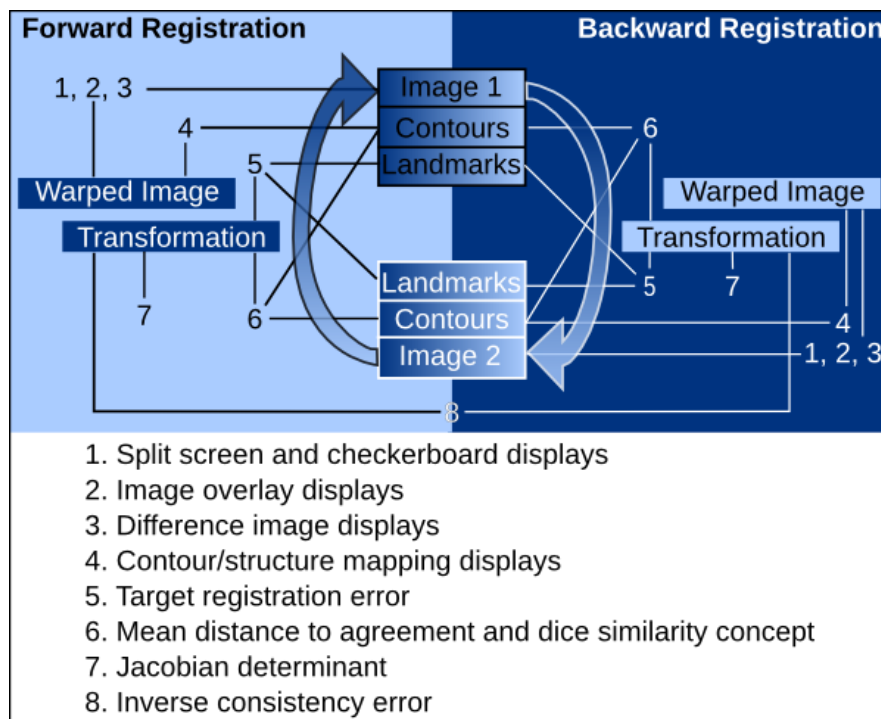


Figure 1 - Schematic presentation of registration QA measures and the required inputs. I.e. for forward measure 4, contours from fixed image and forward warped image are required.

Module Overview

The RegQA module is divided into three parts, as shown in **Error! Reference source not found.**. In the first part user can mark all the inputs. The possible inputs are shown in Figure 1. Additionally, a region of interest (ROI) can be created around the selected segment. The segment has to be first chosen as an input. The ROI is later used in measures 3, 7 and 8. Among two images we can have two registrations (forward and

Slicer RegistrationQA Introduction

backward) and corresponding verifications. User can always switch between the two registrations, which will influence the measure verification execution.

The second part of the module is where different measures can be verified. Buttons for measures become available as soon as necessary inputs are selected in the first part. The overlay of the images has two additional options – inverse color and match window/level. The inverse color will set the overlay images in contrast colors (red and cyan) – all perfectly aligned structures will be therefore displayed in gray. Any deviations will be shown in either red or cyan for fixed or warped (moving) image, respectively. The match window/level checkbox will copy window and level values from the fixed image to the overlaying one. Inverse color and match window/level options are useful when inspecting intra-modality registration.

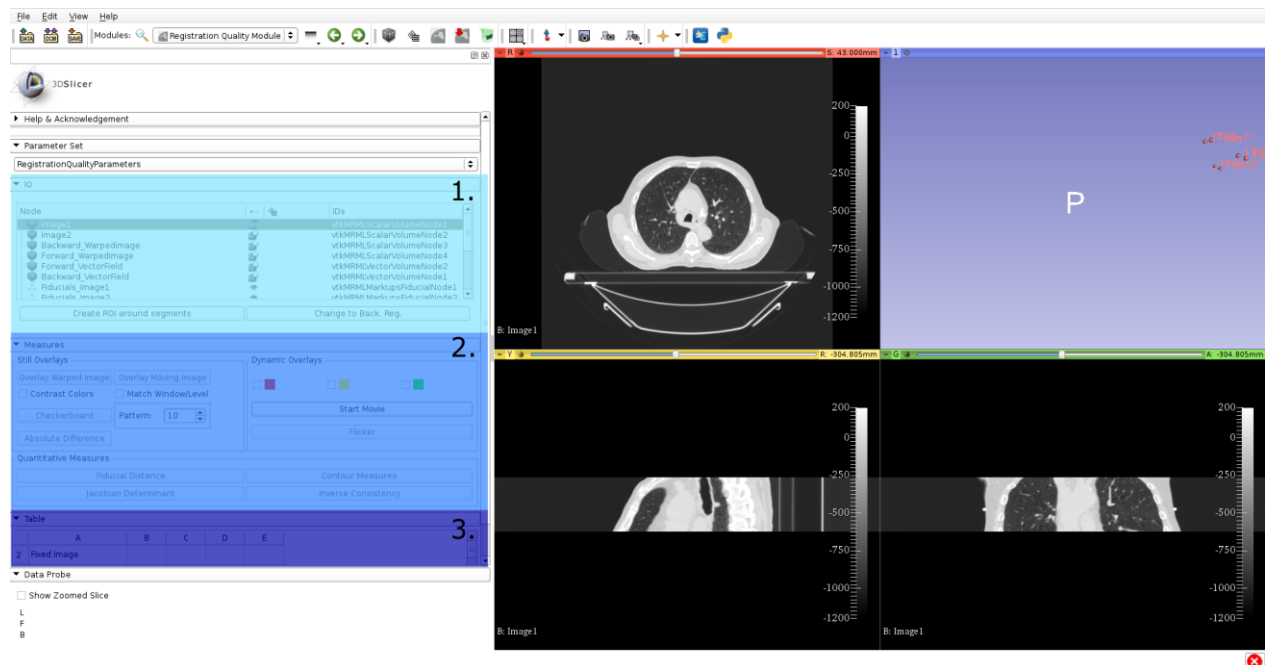


Figure 2 - RegQA module. Three different parts of module are highlighted – input (1), measures (2) and table (3).

The checkerboard and absolute difference correspond to measures 1 and 3. Start movie button automatically increases the slices on the checked views above the button (red, yellow or green). Flicker rapidly switches between foreground and background image, which have to be set previously, either with overlay button or manually.

Slicer RegistrationQA Introduction

Buttons fiducial distance, contour measures, Jacobian determinant and inverse consistency are all quantitative and correspond to measures 5, 6, 7 and 8, respectively. A short overview for them:

1. Fiducial distance calculates distance between fiducials before and after registration. Fiducials (called markups in Slicer) have to be marked in fixed and moving image in the same order – Markup 1 in fixed image has to correspond to Markup 1 in moving image.
2. Contour measures calculate Hausdorff distance and dice metric between contour in fixed and moving image. Currently only one contour (called segment in Slicer) can be selected. It has to be selected in both, fixed and moving image.
3. Jacobian determinant calculates Jacobian determinant on the vector field. Forward and backward Jacobian determinant can be calculated on forward and backward vector fields.
4. Inverse consistency calculates consistency between forward and backward vector field. It has only one value.

The third part on Figure 2 **Error! Reference source not found.** is a table with all the information about the current registration QA – all the inputs with their corresponding nodes and quantitative results from corresponding measures, if they were calculated. The table can also be exported using *Save Data* button in Slicer. An example of the table is shown in Table 1.

Tutorial

Download the data from [here](#)¹ and unzip it in a directory. The images in Image directory represent the two extreme breathing phases. Image1 has been chosen as the reference image and a registration in forward and backward direction has been performed (using Plastimatch module in SlicerRT). The resulting warped images and vector fields are denoted as _warped and _vf, respectively and are stored in RegFiles directory. There are contours and fiducials in Contours and Fiducials directory, respectively, corresponding to both images.

Assigning roles

Load images, warped images, vector fields, contours and fiducials in Slicer. Navigate to *Registration Quality Assurance* module. All inputs should be seen in the input section of module. Next, we have to assign all inputs to their corresponding roles in registration QA, as shown in Figure 1. We'll assign Image 1 to reference image first. Right click on Image 1,

¹ Data downloaded and modified from <https://www.aapm.org/pubs/MPPG/TPS/>

Slicer RegistrationQA Introduction

and then select *Registration QA* -> *Assign node to:* -> *Image* -> *Reference Image*, see Figure 3. If you scroll down in the module, the Image 1 should be seen in the Reference Image row. Repeat the procedure for all other input, assigning them to the right roles. For contours you have to select only GTV63 for both images. If you're unsure of the right nodes, you can check the Table 1 at the end of the document.

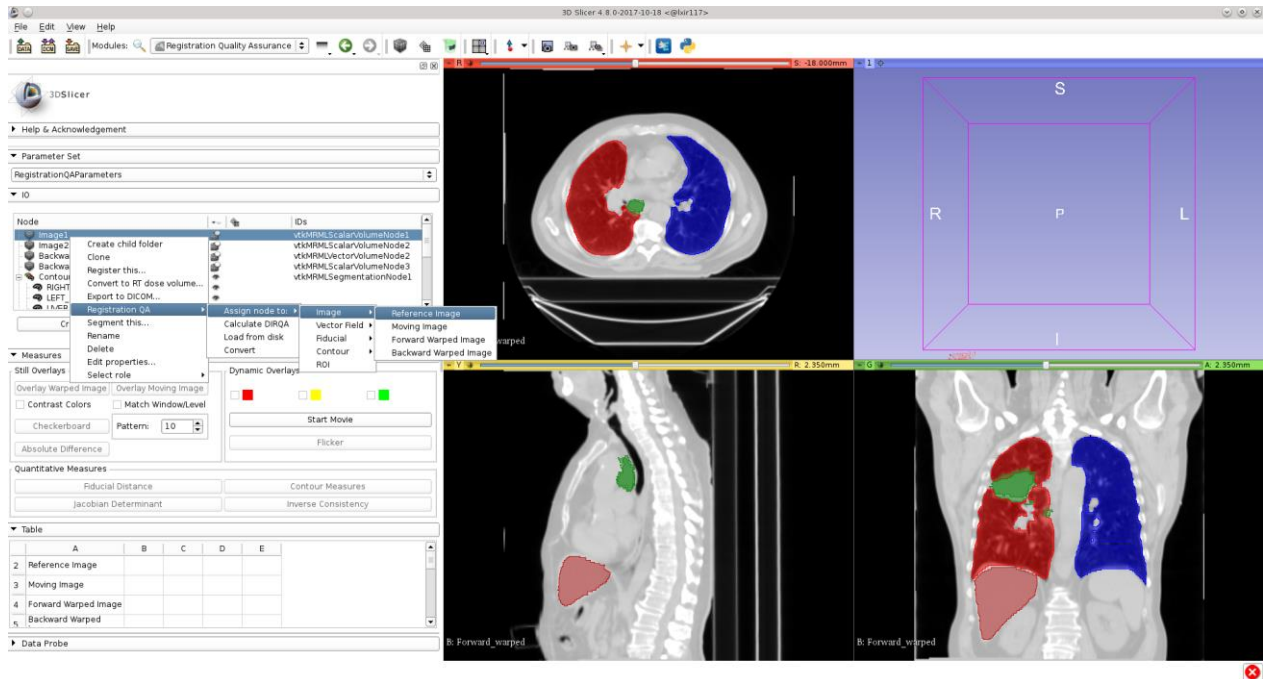


Figure 3 - Assigning roles to specific nodes.

Slicer RegistrationQA Introduction

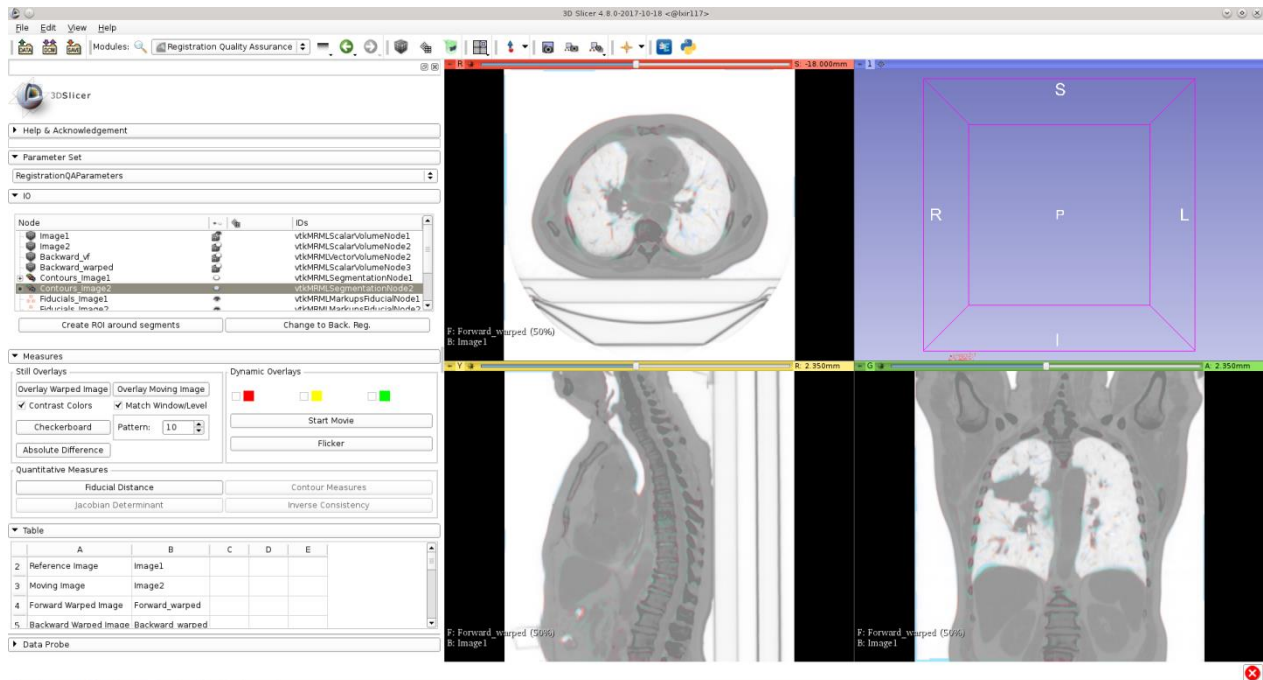


Figure 4 - Overlay with Contrast colors and matched windows and level.

Registration QA checks

Uncheck the eyes next to the `Contours_Image1` and `Contours_Image2`, seen we want to hide their display momentarily. Next, click *Overlay Warped Image* button. The forward warped image will be displayed on top of reference image (Image 1). This corresponds to Measure 2 from Figure 1. Since we are overlaying two CT the difference is not clear. Let's check the *Contrast Colors* and *Match Window/Level* checkboxes. Now click *Overlay Warped Image* button again (Figure 4). Forward warped and reference image are now displayed in opposite colors (cyan and red). Everything that aligns is now shown in gray and deviations are shown in the respective colors. If you click *Overlay Moving Image* button you can see the difference from before registration.

Next we can click *Checkerboard* button. The resulting image shows alternating tiles of forward warped and reference image. This corresponds to Measure 1 and is designed to spot errors in edges.

Absolute difference button will produce an image where each voxel is an absolute difference between values in forward warped and reference image. It corresponds to Measure 3 and mean square error metric, which is commonly used in registration. Details shown in Figure 5.

Slicer RegistrationQA Introduction

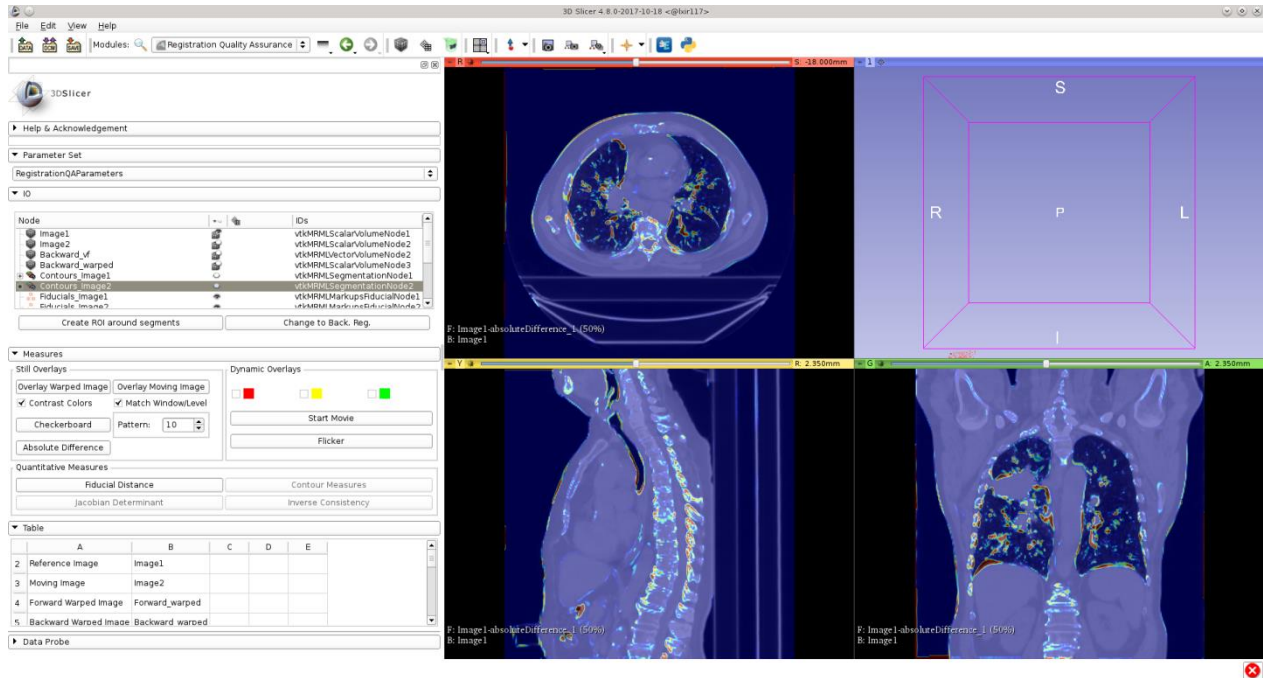


Figure 5 - Absolute difference overlay on top of reference image.

Start Movie and *Flicker* button help user to focus on images. You can select checkbox in front of the red box and click *Start Movie*. The red view should start scrolling slices. You can check other two as well. When you're finished, press *Stop Movie*. Now press *Flicker*. The views will rapidly change between forward warped and reference image. Press *Flicker* again to stop.

To inspect Measure 4 you just have to display only forward warped image and contours from image 1 (reference image). Now you can investigate if contours overlay the desired structures. This usually requires an expert, since contour delineation is not a trivial task.

Next, click *Fiducial Distance*, which corresponds to Measure 5. This calculates distance between fiducials in reference and moving image before and after registration. Results can be found in the table in the bottom part of the module.

Similarly, *Contour Measures* will calculate Hausdorff distance and dice metric between contours before and after registration. You can click the button and investigate results in the table. It corresponds to Measure 6.

Slicer RegistrationQA Introduction

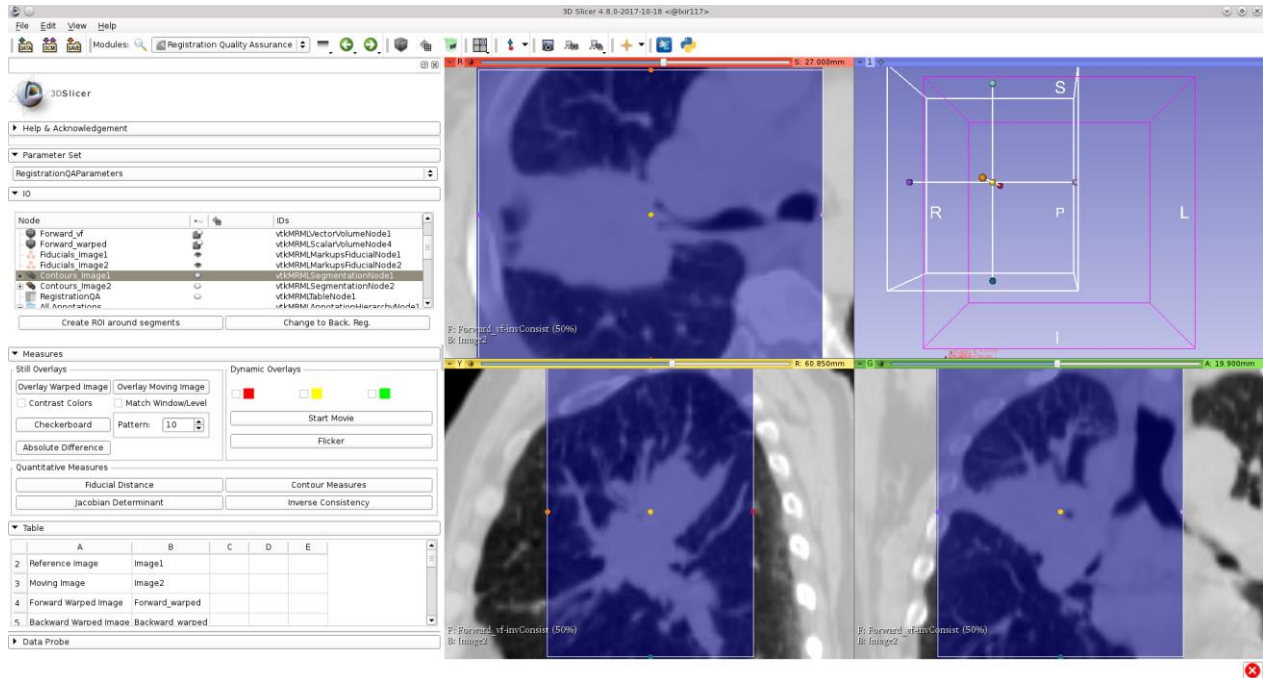


Figure 6 - Inverse consistency calculated inside the ROI around GTV contour.

Jacobian determinant and Inverse Consistency (Measures 7 and 8) will be calculated on the whole CT, if nothing else is specified. This doesn't bring the desired information of registration quality, since a large part of the CT is air. Therefore we want to create a region of interest (ROI) where Jacobian determinant and inverse consistency will be calculated. We can create ROI manually in Slicer and then assign it the right role. Or we can create it automatically by clicking the *Create ROI around segments*, which will create ROI around the contours we assigned. Click this button now and a ROI should appear. Next, click on *Jacobian Determinant* and *Inverse Consistency* buttons. Two images will be displayed over the reference image, one with Jacobian determinant and the other with Inverse Consistency values. The later are normalized to maximum voxel spacing, since inverse consistency error should be smaller than maximum voxel spacing. This means that all inverse consistency values should be below 1. In our case the values are indeed below one, indicating a good registration quality.

Change registration direction

All of the checks explained so far have been for forward registration (except for inverse consistency, which applies to both). We can simply switch which registration direction we want to inspect by clicking *Change to Back. Reg.* Click it now. The button will change to *Change to Forw. Reg.* This enables user to switch between two directions. Repeat now all checks for backward registration.

The table should now be completely filled and be the same as Table 1. Beside in the bottom part of the module you can inspect it in detail in Tables module in Slicer. The table can also be saved to the hard drive and can serve as a registration QA documentation. Alongside should also be screenshots for different measures, similarly as on Figure 7.

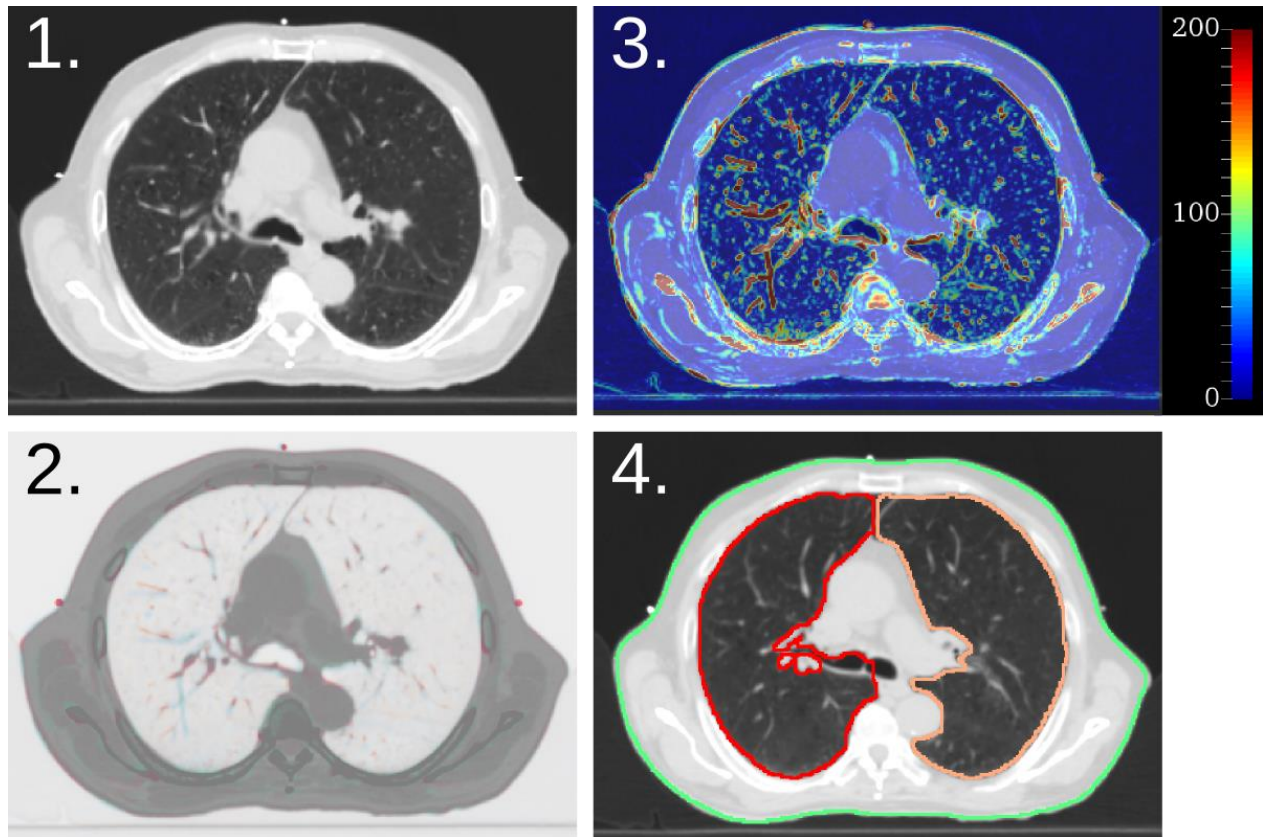


Figure 7 - Examples of Measures 1-4. A checkerboard image with alternating tiles between fixed and warped image is shown on image 1. The image 2 shows overlay of fixed and warped image, where both images are display in contrast colors. An absolute difference between fixed and warped image is shown on image 3. Contour overlay on warped image is shown on image 4.

Create Registration Hierarchy

All the input can be manually loaded; however this takes some time and is tedious work when multiple there's multiple registrations. An auxiliary module called *Create Registration Hierarchy* can create a hierarchy in Slicer, which can be then used to automate the process. The module is written for a specific nomenclature, so you should modify it to serve your needs.

If you open *Create Registration Hierarchy* module and select *Example* from the drop box all the lines should fill up. Replace the text so that it will point to the location on your hard drive where you extracted zip file from previous tutorial. Next click the *Create Subject Hierarchy* button. If you navigate now back to *Registration Quality Assurance*

Slicer RegistrationQA Introduction

module you should see a Subject hierarchy with the name *Example*. Make sure that node under *Parameter Set* is called *ExamplereregQANode*. If not, switch to it. Next right click on the *Example* in hierarchy and click *Registration QA ->Load from disk*. This will load all the files from disk directly into Slicer and assign them the right nodes. With that you can immediately start with specific checks (except for loading of contours which is yet to be implemented).

Additionally, you can right click on *VectorFields* folder in hierarchy view and select *Registration QA ->Calculate Registration QA*, which will perform Measures 7 and 8 on the vector fields. There can also be several vector fields there, such as in 4DCT registration, significantly improving the process. The results for Measures 7 and 8 are also written in two separate tables.

Slicer RegistrationQA Introduction

Table 1 - Input and results from registration QA in from of a table.

Column 1	Column 2	Column 3	Column 4	Column 5
Reference Image	Image1			
Moving Image	Image2			
Forward Warped Image	Forward_warped			
Backward Warped Image	Backward_warped			
Forward vector field	Forward_vf			
Backward vector field	Backward_vf			
Fixed Contour	GTV63			
Moving Contour	GTV63			
Fixed Fiducial	Fiducials_Image1			
Moving Fiducial	Fiducials_Image2			
ROI	GTV63			
Measure:	Max	Min	Mean	STD
AbsDiff (for)	875	0	57.9797	83.7474
AbsDiff (back)	894	0	49.8978	75.8831
Jacobian (for)	1.06519	0.841157	0.940259	0.0435678
Jacobian (back)	1.18093	0.932946	1.06183	0.0501429
InvConsist	0.571222	0.00460278	0.239684	0.103036
	Mean	Haus.	Mean	Haus.
	Before	After	Dice Coeff. Before	Dice Coeff. After
Contour (for)	2.19623	1.27664	0.839715	0.901302
Contour (back)	2.19519	1.13237	0.839845	0.91102
Fiducial	Distance Before	Distance After		
Image1-1	3.80245	2.23676		
Image1-2	4.16145	1.96649		
Image2-1	3.80245	2.15441		
Image2-2	4.16145	1.90811		

References

[1] Brock KK, et al: Use of image registration and fusion algorithms and techniques in radiotherapy: Report of the AAPM Radiation Therapy Committee Task Group No. 132. Medical Physics 44, 2017